



The aSpect Device Engine Python API

aSpect Systems

Oct 30, 2024

CONTENTS

0.1	Getting startet with the Aspect Device Engine Python API	1
0.2	Getting startet with the Aspect Device Engine Python API 2	6
0.3	Getting startet with the Aspect Device Engine Python API 3	8
0.4	Getting startet with the Aspect Device Engine Python API 4	13
0.5	Getting startet with the Aspect Device Engine Python API 5	20
0.6	Getting startet with the Aspect Device Engine Python API 6	22
0.7	Getting startet with the Aspect Device Engine Python API 7	30
0.8	Getting startet with the Aspect Device Engine Python API 8a	30
0.9	Getting startet with the Aspect Device Engine Python API 8b	34

Python Software API Manual



Note

Preliminary Version - This is WIP

Meanwhile visit: [aSpect Systems](#)

[Open Local API documentation](#)

0.1 Getting startet with the Aspect Device Engine Python API

0.1.1 Basics

This is the introductory overview of programming the Aspect Device Engine Python API.

This document is available as pdf and interactive jupyter notebook. The introduction includes the following objectives:

- The few steps to initializing software and hardware
- How to retrieve information about the hardware
- Providing a basic understanding of the structure of the API
- The programming of essential channel parameters
- Setting a voltage and measuring voltage and current

The following spoiler shows a python code snippet and few lines of code that are necessary with the API to generate and measure a voltage:

```
from aspectdeviceengine.enginecore import IdSmuServiceRunner
from aspectdeviceengine.enginecore import IdSmuService, IdSmuBoardModel

# 3 lines of code for the setup
srunner = IdSmuServiceRunner()
mbX1 : IdSmuBoardModel = srunner.get_idsmu_service().get_first_board()
channel1 = mbX1.idsmu2Modules['M1.S1'].smu.channels[1]

# 3 lines of code for configuration and measurement
channel1.enabled = True
channel1.voltage = 2
print(channel1.voltage) # output : ~2.0
```

At the end of this document, this code and some of the background to it should be understandable.

Python imports

There are only a few python imports needed for this introduction. Everthing is imported from the *aspectdeviceengine.enginecore* module. Actually, only the *IdSmuServiceRunner* would be needed since this is the only class that will be instantiated. The other classes are only imported for type hinting. The objects of these types are intantiated by the API services.

```
from aspectdeviceengine.enginecore import IdSmuService, IdSmuServiceRunner, \
↳ IdSmuBoardModel
```

Starting the services and hardware initialization

IdSmuServiceRunner

The *IdSmuServiceRunner* holds the references to the background services. If it goes out of scope, all services are shut down (cleanup processes). The lifetime should therefore be guranteed until the end of the session:

```
srunner = IdSmuServiceRunner()
```

IdSmuService

idSmu devices are detected by the **IdSmuService**. If the *get_first_board()* method is called prior to the detection method, the detection and initialization is performed automatically. This is useful for situations where no specific configuration needs to be done before initialization.

Important note: At the end of a session (be it a jupyter notebook or a python script) the services must be shut down. In the case of the termination of a python script, this happens automatically. When moving from one notebook tutorial to the next, either the kernel must be terminated or the shutdown() method must be executed manually (see last cell in the notebook).

The IdSmuBoardModel

The *IdSmuBoardModel* is the host and (multiside-)controller for idSmu devices.

```
mbX1 : IdSmuBoardModel = srunner.get_idsmu_service().get_first_board()
```

With this single line of code the Hardware is detected and initialized!

Let's print some basic information about the detected devices for this board.

The most relevant information is the **DeviceId** and the device type. The *DeviceId* is used as resource identifier/ locator for the different parts of the hardware. The format is \Mx.Sy.Cz, where x is the mainboard address, y is the device/slot number and z is the channel number.

In the case of the API, the terms *Resource-Id* and *Address* are sysnonyms for the same thing.

```
print(mbX1.device_information)
```

```
+-----+-----+-----+-----+-----+
| DeviceId | IdSmu-Type | Name   | Firmware | Initialized |
+-----+-----+-----+-----+-----+
| M1.S1    | IdSmu2     | M1.S1  | 0x08191f | true        |
+-----+-----+-----+-----+-----+
```

Programming the hardware with the API

idSmu Modules/Devices and Channels

With this system of hierarchical resource localization, each resource is uniquely identifiable, even in a multi-board setup. As we can see (in a single idSmu board setup), there is exactly one device with the address “M1.S1”. An idSmu device (sometimes called module) can contain one or more channels. We can obtain more information about these channels, for example their IDs/addresses:

```
for idSmu2 in mbX1.idSmu2Modules.as_list():
    print(f'idSmu2 ID:      {idSmu2.hardware_id}')
    print(f'idSmu2 name:    {idSmu2.name}')
    print(f'channel IDs:    {idSmu2.channel_ids}')
```

```
idSmu2 ID:      M1.S1
idSmu2 name:    M1.S1
channel IDs:    ['M1.S1.C1', 'M1.S1.C2', 'M1.S1.C3', 'M1.S1.C4']
```

IdSmuDeviceModel

The idSmuModules classes are proxy classes that implement the [] operator for quick access to a device/module of type IdSmuDeviceModel and the as_list() method to get all devices of the same type on the board. There are implementations for the all types of idSmu.

To access a module we can simply use idSmu2Modules['address or name of module']

```
idSmu2 = mbX1.idSmu2Modules['M1.S1']
print(f"The module's id is {idSmu2.hardware_id} and the name is {idSmu2.name}")
```

```
The module's id is M1.S1 and the name is M1.S1
```

Alias names

Devices can be renamed, either programmatically or by so-called parameter settings. Parameter settings are applied during initialization and the changed name can thus be used immediately. (programming via parameter settings is an advanced topic and will not be dealt with here).

The advantage of renaming resources is that you can use an alias name for addressing instead of the rather abstract resource IDs / addresses:

```
idSmu2.name = 'MyFavoriteModule'
# now we can reference the module with the new name
my_fav_module = mbX1.idSmu2Modules['MyFavoriteModule']

print(f"The module's id is still {my_fav_module.hardware_id} and the new name is {my_
↳ fav_module.name}")
```

```
The module's id is still M1.S1 and the new name is MyFavoriteModule
```

The board's device information now lists the new name:

```
print(mbX1.device_information)
```

```
+-----+-----+-----+-----+-----+
| DeviceId | IdSmu-Type | Name           | Firmware | Initialized |
+-----+-----+-----+-----+-----+
| M1.S1    | IdSmu2     | MyFavoriteModule | 0x08191f | true        |
+-----+-----+-----+-----+-----+
```

Descending further into the model hierarchy

The `idSmu-Hardware` is not just a source measurement unit, but combines various hardware components. For example, digital signals can be generated or a RAM memory can be used, depending on the hardware/software support. The software therefore models the hardware in dedicated units. One of the most important elements of the `IdSmuDeviceModel` mentioned above is the unit with which currents and voltages can be generated and measured. These units are called SMU or DPS, depending on the device type.

```
print(idSmu2.smu)
```

```
<aspectdeviceengine.enginecore.IdSmu2DeviceModel.Smu object at 0x000002056E17D9B0>
```

The channel models

The `smu/dps` subcomponents are again proxy objects and contain a `Channels` object.

This `Channels` object implements the `[]` operator for fast access to the channels of the source measurement unit. There is also a `as_list()` method again to iterate over the channels. The objects returned by the channels object are of type `AnalogChannelModel`. This class contains a large part of the methods and properties that you have to deal with in your daily work with the API.

```
for i, channel in enumerate(idSmu2.smu.channels.as_list()):
    print(f"Channel number {i+1} with name {channel.name}"
          f" and identifier {channel.hardware_id}")
    channel.name = f'MyChannel{i+1}'
```

```
Channel number 1 with name M1.S1.C1 and identifier M1.S1.C1
Channel number 2 with name M1.S1.C2 and identifier M1.S1.C2
Channel number 3 with name M1.S1.C3 and identifier M1.S1.C3
Channel number 4 with name M1.S1.C4 and identifier M1.S1.C4
```

After renaming a channel, there are 3 ways to address it:

- through the channel number (starting from 1)
- through the unique channel identifier
- through the channel name

```
print(idSmu2.smu.channels[1].hardware_id)
print(idSmu2.smu.channels["M1.S1.C1"].hardware_id)
print(idSmu2.smu.channels["MyChannel1"].hardware_id)
```

```
M1.S1.C1
M1.S1.C1
M1.S1.C1
```

Important Note: Attempting to assign the same name to two different channels leads to an exception. Channel names must be unique. The reason is that the engine accepts names as resource identifiers for many

operations. If the user were given the option to overwrite this name, bugs that are difficult to identify would be possible:

```
try:
    idSmu2.smu.channels[2].name="MyChannel1"
except Exception as e:
    print(f"Exception: {str(e)}")
```

```
Exception: The alias name MyChannel1 is already associated with the id M1.S1.C1
```

Parameterization and measurements

Now we have a reference to the channel object and can finally parameterize it and take measurements. A channel must be active so that a voltage (or a current) can be output or meaningful measurements can be made.

Let's check if the channel is enabled, and if it is not, enable it:

```
channel1 = idSmu2.smu.channels["MyChannel1"]
print(f'Channel enabled? {channel1.enabled}')
if not channel1.enabled:
    channel1.enabled = True
print(f'Channel enabled? {channel1.enabled}')
```

```
Channel enabled? False
Channel enabled? True
```

Measuring voltage and current

The quickest and easiest way to measure a voltage or a current are the properties voltage and current

```
print(f'Measured voltage: {channel1.voltage:6f}')
print(f'Measured current: {channel1.current:6f}')
```

```
Measured voltage: 0.000000
Measured current: 0.000002
```

Setting voltage and current

The quickest and easiest way to set a voltage or a current are the setters voltage and current

```
channel1.voltage = 3.14
print(f'Measured voltage: {channel1.voltage:6f}')

# Setting a current is only usefull if there is a load at the outputs
# wheras voltages can be measured on an open output
# channel1.current = 1E-3
# print(f'Measured current: {channel1.current:6f}')
```

```
Measured voltage: 3.141083
```

Bonus: Getting the maximum output voltage and output current ratings

```
vMin, vMax, iMin, iMax = channel1.output_ranges
print(f'Output voltage range: [{vMin:6f}, {vMax:6f}] V')
print(f'Output current range: [{iMin:6f}, {iMax:6f}] A')
```

```
Output voltage range: [-11.000000, 10.999664] V
Output current range: [-0.075000, 0.074998] A
```

```
srunner.shutdown()
```

0.2 Getting started with the Aspect Device Engine Python API 2

0.2.1 Multiple Measurements

This is the second introductory overview of programming the Aspect Device Engine Python API.

This document is available as pdf and interactive jupyter notebook. The introduction includes the following objectives:

- Get familiar with multiple measurements at multiple channels

```
from aspectdeviceengine.enginecore import IdSmuService, IdSmuServiceRunner, \
↳ IdSmuBoardModel
import plotly.graph_objects as go
import numpy as np
srunner = IdSmuServiceRunner()
```

```
mbX1 : IdSmuBoardModel = srunner.get_idsmu_service().get_first_board()
```

```
channel1 = mbX1.idSmu2Modules['M1.S1'].smu.channels["M1.S1.C1"]
channel2 = mbX1.idSmu2Modules['M1.S1'].smu.channels["M1.S1.C2"]
channel3 = mbX1.idSmu2Modules['M1.S1'].smu.channels["M1.S1.C3"]
channel4 = mbX1.idSmu2Modules['M1.S1'].smu.channels["M1.S1.C4"]
print(channel1)
```

```
<aspectdeviceengine.enginecore.AD5522ChannelModel object at 0x00000290455AB8B0>
```

0.2.2 Repeated measurements

While measurements can be repeated in the software, e.g. in a loop, as shown in the previous sections, the measurements can also be repeated directly on the hardware.

This has the advantage that the code is much more compact and that the measurements run much faster.

The **measure_voltages()** method takes one parameter: the number of repetitions for the measurement.

```
channel1.enabled = True
channel1.voltage = 1.14
number_of_measurements = 50
voltages_ch1 = channel1.measure_voltages(number_of_measurements)
print(voltages_ch1[1:10])
```

```
[1.138671875, 1.138671875, 1.138336181640625, 1.138336181640625, 1.138336181640625,
↪ 1.138336181640625, 1.137664794921875, 1.138336181640625, 1.137664794921875]
```


Repeat executing the next cell (Ctrl+Enter) and you will see how the below plot updates.

The plot stays close to 1.4V. Only the usual very small fluctuations due to noise can be seen.

```
number_of_measurements = 50
voltages_ch1 = channel1.measure_voltages(number_of_measurements)

x_ = np.linspace(0, number_of_measurements, number_of_measurements)
# Create traces
fig = go.Figure()
fig.add_trace(go.Scatter(x=x_, y=voltages_ch1, mode='lines+markers', name='voltages_
↪ch1'))

fig.update_layout( title={'text': f"Single channel measurement, \
std.dev: {np.std(voltages_ch1):.6f}, mean: {np.mean(voltages_ch1):.6f}", \
'y':0.9, 'x':0.5, 'xanchor': 'center', 'yanchor': 'top'}, \
axis_title='Measurement repetition number', yaxis_title='Voltage [V]', \
margin=dict(l=20, r=20, t=55, b=20))
fig
#uncomment in pure python script:
#fig.show()
```

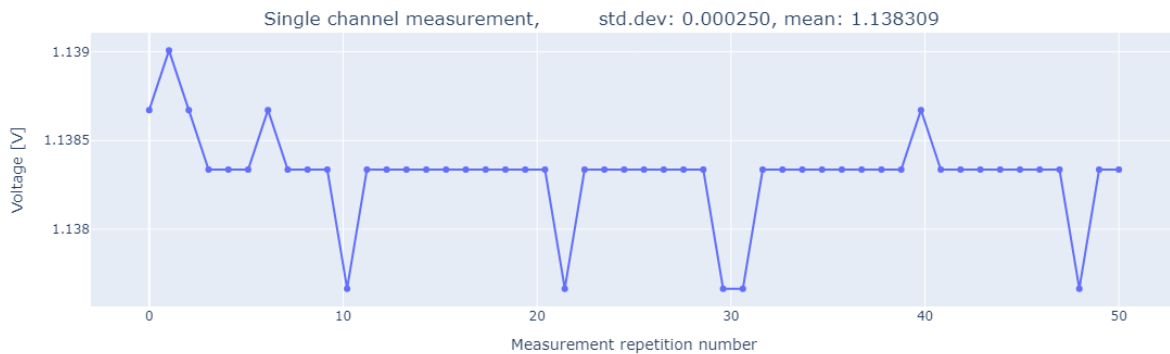

Now we activate the other channels with an offset of one volt each.

```
channel2.enabled = True
channel2.voltage = 2.14
channel3.enabled = True
channel3.voltage = 3.14
channel4.enabled = True
channel4.voltage = 4.14
```

The measurements for each channel are run sequentially and the results are displayed in a plot:

```
voltages_ch1 = channel1.measure_voltages(50)
voltages_ch2 = channel2.measure_voltages(50)
voltages_ch3 = channel3.measure_voltages(50)
voltages_ch4 = channel4.measure_voltages(50)

x_ = np.linspace(0, 50, 50)
# Create traces
fig = go.Figure()
```

(continues on next page)

(continued from previous page)

```

fig.add_trace(go.Scatter(x=x_, y=voltages_ch1,
                        mode='lines+markers',
                        name='voltages_ch1'))
fig.add_trace(go.Scatter(x=x_, y=voltages_ch2,
                        mode='lines+markers',
                        name='voltages_ch2'))
fig.add_trace(go.Scatter(x=x_, y=voltages_ch3,
                        mode='lines+markers',
                        name='voltages_ch3'))
fig.add_trace(go.Scatter(x=x_, y=voltages_ch4,
                        mode='lines+markers',
                        name='voltages_ch4'))
fig.update_layout( title={'text': "4 channel measurement", 'y':0.9, 'x':0.5,
↪ 'xanchor': 'center', 'yanchor': 'top'},
                  axis_title='Measurement repetition number', yaxis_title='Voltage_
↪ [V]', margin=dict(l=20, r=20, t=55, b=20))
fig
#uncomment in pure python script:
#fig.show()

```



Do not forget to shut down the services before proceeding:

```
srunner.shutdown()
```

0.3 Getting started with the Aspect Device Engine Python API 3

0.3.1 Multi-side and high performance

This is the third introductory overview of programming the Aspect Device Engine Python API.

This document is available as pdf and interactive jupyter notebook. The introduction includes the following objectives:

- Understanding the parallel or multi-side programming of idSMU resources

Introduction

In the previous chapters, the idSMU channels were always programmed individually and sequentially.

This is the simplest and quickest approach to programming the hardware. This approach is sufficient in cases where a high-performance application is not important.

With its multi-threading approach, the aspect device engine also offers the possibility to programme and measure resources efficiently in parallel.

```
from aspectdeviceengine.enginecore import IdSmuService, IdSmuServiceRunner, \
    IdSmuBoardModel, MeasurementMode
import plotly.graph_objects as go
import numpy as np
srunner = IdSmuServiceRunner()
```

```
mbX1 : IdSmuBoardModel = srunner.get_idsmu_service().get_first_board()
mbX1.is_board_initialized()
```

```
True
```

The IdSmuBoardModel class is the host and controller for the detected hardware on an idSMU board.

A board can, for example, contain a single idSMU module (MbX-1) or 16 (MbX-16). Parallelisation is optimised for one board, i.e. all theoretically possible 160 channels of an MbX-16 board could be set to 0V at once and the engine would attempt to carry out this process as efficiently as possible.

Parallelism between several boards is also given in the sense that the programming of a resource is a non-blocking process running in the background.

The most frequently used methods for parameterising and measuring a channel are also available as board methods.

For example, the property `.voltage` of a channel corresponds to the `set_voltages()` method of a board.

The plural in the method name already indicates that several channels can be programmed simultaneously here.

```
# The anatomy of a board method usually consists of the combination parameter-value,
# followed by a list of device or channel names (or identifiers)
# We could set a voltage for each channel separately:
mbX1.set_voltages(3.14, ["M1.S1.C1"])
mbX1.set_voltages(4.13, ["M1.S1.C2"])
# or in parallel:
mbX1.set_voltages(3.14, ["M1.S1.C1", "M1.S1.C2"])

# After assigning some alias names...
mbX1.idSmu2Modules['M1.S1'].smu.channels["M1.S1.C1"].name = "channel1"
mbX1.idSmu2Modules['M1.S1'].smu.channels["M1.S1.C2"].name = "channel2"
# ...the new names can be used in the board methods:
mbX1.set_voltages(3.14, ["channel1", "channel2"])
# Changing the measurement mode and enabling the channels is done in a similar way:
mbX1.set_measurement_modes(MeasurementMode.vsense, ["channel1", "channel2"])
mbX1.set_enable_channels(True, ["channel1", "channel2"] )
```

Parallel measurements in detail

```
smuchannel1 = mbX1.idSmu2Modules['M1.S1'].smu.channels['channel1']
smuchannel2 = mbX1.idSmu2Modules['M1.S1'].smu.channels['channel2']
```

```
print(smuchannel1.enabled, smuchannel2.enabled)
```

```
True True
```

```
print(smuchannel1.voltage, smuchannel2.voltage)
```

```
3.13873291015625 3.14007568359375
```

We have set and measured the output voltage using the simple and intuitive method via channel properties, as presented in the last tutorials. These measurements run sequentially. Next, we will learn about the high-performance parallel method. The `IdSmuBoardModel` has numerous methods for programming and measuring the resources it manages. One of these methods is `set_voltages()`, which can be used to set any number of channels together to the same voltage value.

```
mbX1.set_voltages(voltage=2.5, channel_names=["channel1", "channel1"])
# for SMU-based modules, set_currents() can be called as well
```

We can still use the “channel-method” to query the voltages:

```
print(smuchannel1.voltage, smuchannel2.voltage)
```

```
2.099090576171875 3.14007568359375
```

Synchronous and asynchronous measurements

When performing a measurement via the board controller (*IdSmuBoardModel*) on one or more channels, you have the choice of either waiting for the result until the measurement has been performed or letting the measurement run in the background. The latter is particularly useful for triggered measurements or for GUI applications or other asynchronous tasks.

This is set with the `wait_for_result` parameter of the `measure_channels()` method.

The `measure_channels()` method also offers the advantage of being able to set other parameters, such as the sample count or the number of measurements. Let’s perform a synchronous measurement with 2 repetitions on 2 channels:

```
measresults = mbX1.measure_channels(wait_for_result=True, sample_count=1,
↳ repetitions=2, channel_names=["channel1", "channel2"])
```

Important Note:

The user is responsible to provide a valid list of channel names/identifiers. In addition to the valid names, care must also be taken to ensure that the channel names in the list are unique without duplicate entries. Otherwise errors will occur!

Measurement Results

The result of a measurement via a board is a vector (list) of measurement results.

Each of these measurement results relates to an idSmu module/device.

As the `measure_channels()` method can be used to measure several channels on several modules quasi-parallel, this result list can contain more than one entry. As we only measured on one module on two channels, the list only has one entry:

```
print(len(measresults))
measresult0 = measresults[0]
type(measresult0)
```

```
1
```

```
aspectdeviceengine.enginecore.ReadAdcCommandIdSmuResult
```

The elements of the result are of type `ReadAdcCommandIdSmuResult`. Useful properties are `device_id`, `channel_ids`, `channel_names`, `execution_time`:

```
print(f'The results come from the measurement on the device with the id {measresult0.
↪device_id}')
print(f'The total execution time (including data transfer via usb etc) was
↪{measresult0.execution_time} microseconds.')
```

```
The results come from the measurement on the device with the id M1.S1
The total execution time (including data transfer via usb etc) was 669_
↪microseconds.
```

To further simplify the assignment of the results to the measured resources, we can query the IDs of the channels or their names (the names must have been set before execution, see above)

```
print(measresult0.channel_ids, measresult0.channel_names)
```

```
StringList[M1.S1.C1, M1.S1.C2] StringList[channel1, channel2]
```

There are now several ways to obtain the result for a specific channel (as numpy array), which are all equivalent:

```
print(measresult0.get_float_values("M1.S1.C1"))
print(measresult0["M1.S1.C1"])
print(measresult0["channel1"])
print(measresult0[measresult0.channel_names[0]])

# example of averaging the results for each channel:
for channel_name in measresult0.channel_names:
    print(f'{channel_name} with an average value of {np.mean(measresult0[channel_
↪name]) :.4f} V')
```

```
[2.09875488 2.09942627]
[2.09875488 2.09942627]
[2.09875488 2.09942627]
[2.09875488 2.09942627]
channel1 with an average value of 2.0991 V
channel2 with an average value of 3.1399 V
```

Asynchronous measurements

With asynchronous measurements, there is no waiting for a measurement. The commanded measurements are started in the background in high-performance C++ threads. The measurement results can be retrieved in python at any time. The result of the non-waiting `measure_channels()` method is therefore an empty array as shown below:

```
measresults = mbX1.measure_channels(wait_for_result=False, sample_count=1,
    repetitions=2, channel_names=["channel1", "channel2"])
print(measresults)
```

```
[]
```

The `get_measurement_results_for_channel()` method returns the result of at least the specified channel. If several channels are measured simultaneously on one device, as in this example, all results are returned (a device cannot return the results separately for each channel). Since *channel1* and *channel2* channels are on the same device, an object is returned that bundles the results for these two channels. It is of the same type as the element already from the array after the synchronous call of the measurement. To recognize it, we call it “measresult0” again

```
measresult0 = mbX1.get_measurement_results_for_channel("channel1")
```

```
print(measresult0["channel1"])
print(measresult0.timecode)
```

```
[2.10043335 2.10043335]
[0 0]
```

Timecode generation during measurements

In addition to the actual measurement data, a time code can be generated that tracks the exact time at which a measurement was started. The following method is used to activate this (after restarting a device, the default is *disabled*):

```
mbX1.enable_timecode("M1.S1")
```

The timecode is counted in multiples of 10ns (100Mhz clock). It is generated by counter that always runs when enabled, not only when commands are sent. We subtract the first value from the array for an offset of zero and divide by 100 to get the time in units of microseconds.

```
measresults = mbX1.measure_channels(wait_for_result=False, sample_count=1,
    repetitions=20, channel_names=["channel1", "channel2"])
measresult0 = mbX1.get_measurement_results_for_channel("channel1")

timecode = (measresult0.timecode-measresult0.timecode[0])/100
```

Now we can display the measurement results in a plot over an axis that displays the normalized time vs the measurement results

```
x_ = timecode
# Create traces
fig = go.Figure()
fig.add_trace(go.Scatter(x=x_, y=measresult0["channel1"],
    mode='lines+markers',
    name='ch1'))
fig.add_trace(go.Scatter(x=x_, y=measresult0["channel2"],
    mode='lines+markers',
```


(continues on next page)

(continued from previous page)

```

        name='ch2'))
fig.update_layout( title={'text': "Parallel 2 channel measurement", 'y':0.9, 'x':0.
↪5, 'xanchor': 'center', 'yanchor': 'top'},
                  axis_title='Time [us]', yaxis_title='Voltage [V]',
↪margin=dict(l=20, r=20, t=55, b=20))
fig
#uncomment in pure python script:
#fig.show()

```

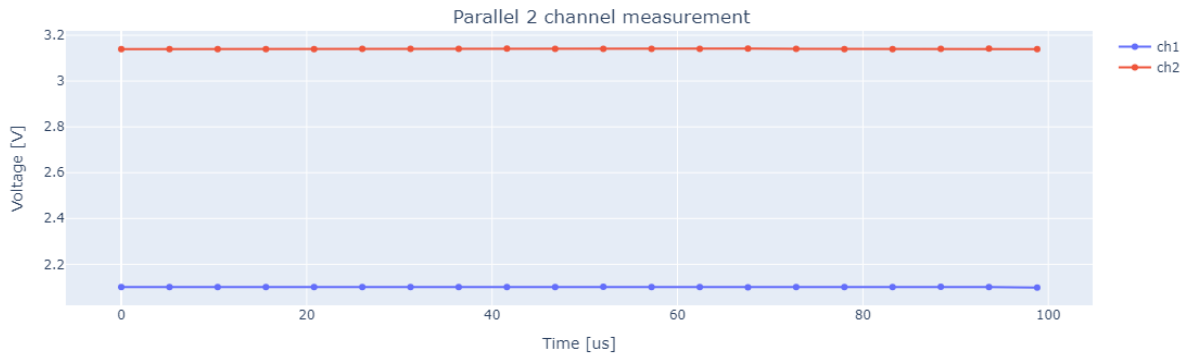

Do not forget to shut down the services before proceeding:

```
srunner.shutdown()
```

0.4 Getting startet with the Aspect Device Engine Python API 4

0.4.1 Current Range and Clamps

This is the forth introductory overview of programming the Aspect Device Engine Python API.

This document is available as pdf and interactive jupyter notebook. The introduction includes the following objectives:

- Understanding the current range and clamps of idSMU devices

```

from aspectdeviceengine.enginecore import IdSmuService, IdSmuServiceRunner,
↪IdSmuBoardModel, MeasurementMode, SmuCurrentRange, DpsCurrentRange, CurrentRange
import plotly.graph_objects as go
from plotly.subplots import make_subplots
import numpy as np
srunner = IdSmuServiceRunner()

```

```

mbX1 : IdSmuBoardModel = srunner.get_idsmu_service().get_first_board()
print(mbX1.is_board_initialized())
smu_channel = mbX1.idSmu2Modules['M1.S1'].smu.channels["M1.S1.C1"]
smu_channel.name = "ch1"
smu_channel.enabled = True

```

```
True
```

Python can always be used to analyse the properties and methods of an object.
We search for current range and clamp in the list, the topics of this tutorial:

```
# Let's list all properties and methods of the smu channel object
[m for m in dir(smu_channel) if not m.startswith('_')]
```

```
['autorange',
 'clamp_enabled',
 'clamp_high_value',
 'clamp_low_value',
 'current',
 'current_range',
 'enabled',
 'hardware_id',
 'measure_current',
 'measure_currents',
 'measure_voltage',
 'measure_voltages',
 'name',
 'output_ranges',
 'set_clamp_high_value',
 'set_clamp_low_value',
 'set_current',
 'set_name',
 'set_voltage',
 'voltage']
```

Current Range

As an API user, you have the option of manually setting different current ranges.

It often makes sense to select a current range that is as close as possible to the current flowing through the DUT.

For small currents, the smallest possible matching current range should be selected.

This increases the accuracy when measuring a current. The two possible device types of an idSMU (SMU and DPS) each have different current ranges.

Setting and querying the current ranges

The enum class for the current ranges only contains the respective valid members at channel level for SMU or DPS.

The enum at board level is a union of both.

In the latter case, the user is responsible for not selecting the wrong value.

```
print(list(SmuCurrentRange.__members__))
print(list(DpsCurrentRange.__members__))
print(list(CurrentRange.__members__))
```

```
['Range_5uA', 'Range_20uA', 'Range_200uA', 'Range_2mA', 'Range_70mA']
['Range_25uA', 'Range_250uA', 'Range_2500uA', 'Range_25mA', 'Range_500mA', 'Range_
↪1200mA']
['Range_5uA', 'Range_20uA_SMU', 'Range_200uA_SMU', 'Range_2mA_SMU', 'Range_70mA_SMU
↪', 'Range_25uA_DPS', 'Range_250uA_DPS', 'Range_2500uA_DPS', 'Range_25mA_DPS',
↪'Range_500mA_DPS', 'Range_1200mA_DPS']
```

The default value (after initialization) of a SMU based device is 70mA (and 500mA for a DPS device). The `current_range` property of a channel or the `get_current_range()` method of a board can be used to query the

current range:

```
print(smu_channel.current_range)
print(mbX1.get_current_range("ch1"))
```

```
SmuCurrentRange.Range_70mA
CurrentRange.Range_70mA_SMU
```

The corresponding setter and board method are `current_range` and `set_current_ranges()`

```
smu_channel.current_range = SmuCurrentRange.Range_2mA
print(mbX1.get_current_range("ch1"))
mbX1.set_current_ranges(CurrentRange.Range_200uA_SMU, ["ch1"])
print(smu_channel.current_range)
```

```
CurrentRange.Range_2mA_SMU
SmuCurrentRange.Range_200uA
```

More about current ranges and autoranging in the next tutorial.

Current and voltage clamps

An SMU or DPS channel has clamps. If a voltage is forced, the clamp is a current clamp.

If a current is forced (SMU types only), a voltage is clamped.

The clamp is active by default. We can check this with the 'clamp_enabled' property:

```
print(f'Is clamp enabled? {smu_channel.clamp_enabled}')
```

```
Is clamp enabled? True
```

In voltage force mode, the default clamp value is 70mA

```
print(f'The default upper clamp value is {smu_channel.clamp_high_value}A')
```

```
The default upper clamp value is 0.07A
```

To see the clamp in action, the voltage is swept across an LED and the current is measured (more on IV sweeps in the next chapters).


The clamp is first set to 100uA and the current range is adjusted to this current range:

```
currents = np.zeros(40); voltages = np.zeros(40)
smu_channel.clamp_enabled = True

smu_channel.clamp_high_value=0.0001
# an alternative to set lower and upper clamp together is to use the board method:
# mbX1.set_clamps_low_and_high_values(-0.001, 0.0001, ["ch1"])
smu_channel.current_range = SmuCurrentRange.Range_200uA
for i in range(1,40):
    smu_channel.voltage = 2.0+0.2*i
    currents[i] = smu_channel.current
    voltages[i] = smu_channel.voltage
```

As can be seen in the plot, the clamp becomes active slightly above the set maximum. (Deviations can occur if the clamp is not calibrated)

```
fig = make_subplots(rows=1, cols=2, subplot_titles=("LED current", "LED voltage"))
fig.add_trace(go.Scatter(x=np.arange(1,40), y=currents,
                        mode='lines+markers+text',
                        textposition="top center",
                        name='LED current'), row=1, col=1)
fig.add_trace(go.Scatter(x=np.arange(1,40), y=voltages,
                        mode='lines+markers+text',
                        textposition="top center",
                        name='LED voltage'), row=1, col=2)
fig.update_layout(height=600)
fig.update_xaxes(title_text="Step", row=1, col=1);fig.update_yaxes(title_text=
↪ "Current [A]", row=1, col=1)
fig.update_xaxes(title_text="Step", row=1, col=1);fig.update_yaxes(title_text=
↪ "Coltage [V]", row=1, col=2)
fig
#uncomment in pure python script:
#fig.show()
```

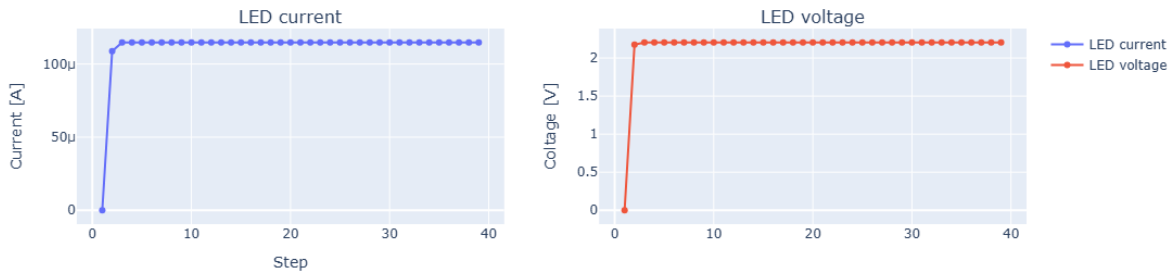

The clamp will now be disabled and the current measurement will be repeated:

```
currents = np.zeros(40);voltages = np.zeros(40)

smu_channel.clamp_high_value=0.0001
smu_channel.clamp_enabled = False
smu_channel.current_range = SmuCurrentRange._2mA
for i in range(1,40):
    smu_channel.voltage = 2.0+0.2*i
    currents[i] = smu_channel.current
    voltages[i] = smu_channel.voltage
```

```
fig = make_subplots(rows=1, cols=2, subplot_titles=("LED current", "LED voltage"))
fig.add_trace(go.Scatter(x=np.arange(1,40), y=currents,
                        mode='lines+markers+text',
                        textposition="top center",
                        name='LED current'), row=1, col=1)
fig.add_trace(go.Scatter(x=np.arange(1,40), y=voltages,
                        mode='lines+markers+text',
                        textposition="top center",
                        name='LED voltage'), row=1, col=2)
fig.update_layout(height=600)
fig.update_xaxes(title_text="Step", row=1, col=1);fig.update_yaxes(title_text=
↪ "Current [A]", row=1, col=1)
fig.update_xaxes(title_text="Step", row=1, col=1);fig.update_yaxes(title_text=
```


(continues on next page)

(continued from previous page)

```

↪ "Coltage [V]", row=1, col=2)
fig
#uncomment in pure python script:
#fig.show()

```

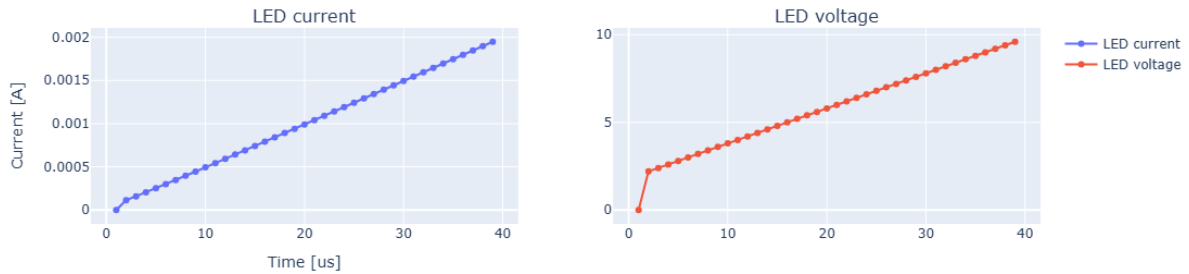

```

srunner.shutdown()

```

```

measresults = mbX1.measure_channels(wait_for_result=True, sample_count=1,
↪ repetitions=2, channel_names=["channel1", "channel2"])

```

Important Note:

The user is responsible to provide a valid list of channel names/identifiers. In addition to the valid names, care must also be taken to ensure that the channel names in the list are unique without duplicate entries. Otherwise errors will occur!

Measurement Results

The result of a measurement via a board is a vector (list) of measurement results.

Each of these measurement results relates to an idSmu module/device.

As the `measure_channels()` method can be used to measure several channels on several modules quasi-parallel, this result list can contain more than one entry. As we only measured on one module on two channels, the list only has one entry:

```

print(len(measresults))
measresult0 = measresults[0]
type(measresult0)

```

```

1

```

```

aspectdeviceengine.enginecore.ReadAdcCommandIdSmuResult

```

The elements of the result are of type `ReadAdcCommandIdSmuResult`. Useful properties are `device_id`, `channel_ids`, `channel_names`, `execution_time`:

```

print(f'The results come from the measurement on the device with the id {measresult0.
↪ device_id}')
print(f'The total execution time (including data transfer via usb etc) was
↪ {measresult0.execution_time} microseconds.')

```

```
The results come from the measurement on the device with the id M1.S1
The total execution time (including data transfer via usb etc) was 669_
↳microseconds.
```

To further simplify the assignment of the results to the measured resources, we can query the IDs of the channels or their names (the names must have been set before execution, see above)

```
print(m measresult0.channel_ids, measresult0.channel_names)
```

```
StringList[M1.S1.C1, M1.S1.C2] StringList[channel1, channel2]
```

There are now several ways to obtain the result for a specific channel (as numpy array), which are all equivalent:

```
print(m measresult0.get_float_values("M1.S1.C1"))
print(m measresult0["M1.S1.C1"])
print(m measresult0["channel1"])
print(m measresult0[measresult0.channel_names[0]])

# example of averaging the results for each channel:
for channel_name in measresult0.channel_names:
    print(f'{channel_name} with an average value of {np.mean(m measresult0[channel_
↳name]) :.4f} V')
```

```
[2.09875488 2.09942627]
[2.09875488 2.09942627]
[2.09875488 2.09942627]
[2.09875488 2.09942627]
channel1 with an average value of 2.0991 V
channel2 with an average value of 3.1399 V
```

With asynchronous measurements, there is no waiting for a measurement. The commanded measurements are started in the background in high-performance C++ threads. The measurement results can be retrieved in python at any time. The result of the non-waiting `measure_channels()` method is therefore an empty array as shown below:

```
measresults = mbX1.measure_channels(wait_for_result=False, sample_count=1,
↳repetitions=2, channel_names=["channel1", "channel2"])
print(m measresults)
```

```
[]
```

The `get_measurement_results_for_channel()` method returns the result of at least the specified channel. If several channels are measured simultaneously on one device, as in this example, all results are returned (a device cannot return the results separately for each channel). Since *channel1* and *channel2* channels are on the same device, an object is returned that bundles the results for these two channels. It is of the same type as the element already from the array after the synchronous call of the measurement. To recognize it, we call it “measresult0” again

```
measresult0 = mbX1.get_measurement_results_for_channel("channel1")
```

```
print(m measresult0["channel1"])
print(m measresult0.timecode)
```

```
[2.10043335 2.10043335]
[0 0]
```

Timecode generation during measurements

In addition to the actual measurement data, a time code can be generated that tracks the exact time at which a measurement was started. The following method is used to activate this (after restarting a device, the default is *disabled*):

```
mbX1.enable_timecode("M1.S1")
```


The timecode is counted in multiples of 10ns (100Mhz clock). It is generated by counter that always runs when enabled, not only when commands are sent. We subtract the first value from the array for an offset of zero and divide by 100 to get the time in units of microseconds.

```
measresults = mbX1.measure_channels(wait_for_result=False, sample_count=1,
    repetitions=20, channel_names=["channel1", "channel2"])
measresult0 = mbX1.get_measurement_results_for_channel("channel1")

timecode = (measresult0.timecode-measresult0.timecode[0])/100
```

Now we can display the measurement results in a plot over an axis that displays the normalized time vs the measurement results

```
x_ = timecode
# Create traces
fig = go.Figure()
fig.add_trace(go.Scatter(x=x_, y=measresult0["channel1"],
    mode='lines+markers',
    name='ch1'))
fig.add_trace(go.Scatter(x=x_, y=measresult0["channel2"],
    mode='lines+markers',
    name='ch2'))
fig.update_layout(title={'text': "Parallel 2 channel measurement", 'y':0.9, 'x':0.
    5, 'xanchor': 'center', 'yanchor': 'top'},
    axis_title='Time [us]', yaxis_title='Voltage [V]',
    margin=dict(l=20, r=20, t=55, b=20))
fig
#uncomment in pure python script:
#fig.show()
```

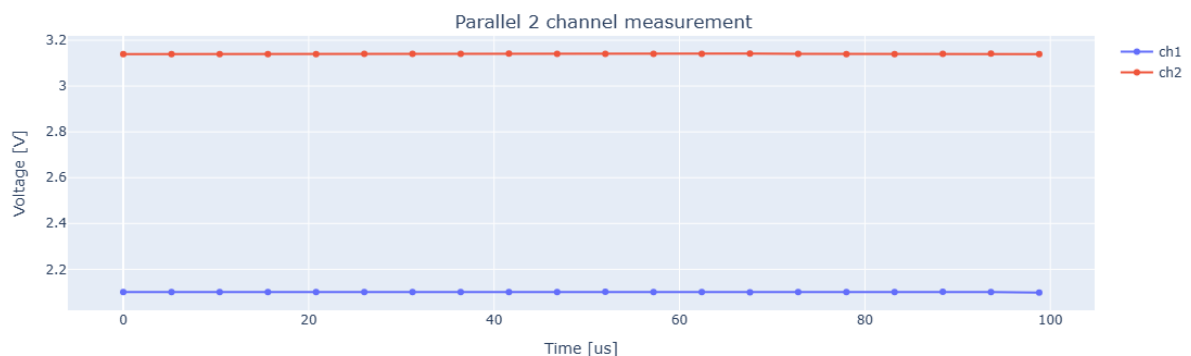

Do not forget to shut down the services before proceeding:

```
srunner.shutdown()
```

0.5 Getting started with the Aspect Device Engine Python API 5

0.5.1 Simple software IV sweeps and autoranging

This is the fifth introductory overview of programming the Aspect Device Engine Python API.

This document is available as pdf and interactive jupyter notebook. The introduction includes the following objectives:

- Performing IV sweeps
- Working with autoranging in IV-sweeps

Please note: This chapter shows how to perform simple IV sweeps in software.

For more advanced IV sweeps, see the next chapter 'List Sweeps'.

```
from aspectdeviceengine.enginecore import IdSmuService, IdSmuServiceRunner, \
↳ IdSmuBoardModel, MeasurementMode, FunctionGeneratorType, generate_function_
↳ generator_data, CurrentRange
from aspectdeviceengine.enginecore import IdSmuBoardModel, \
↳ ListSweepChannelConfiguration, ListSweep
import plotly.graph_objects as go
import numpy as np
srunner = IdSmuServiceRunner()
```

A few preparatory measures, nothing new here:

```
# Initialization and preparation:
mbX1 : IdSmuBoardModel = srunner.get_idsmu_service().get_first_board()
print(mbX1.is_board_initialized())
device_id = "M1.S1"
channel1_id = f'{device_id}.C1'
channel2_id = f'{device_id}.C3'
channel_list = [channel1_id, channel2_id]

mbX1.set_voltages(1, channel_list)
mbX1.set_enable_channels(True, channel_list )
mbX1.set_measurement_modes(MeasurementMode.isense, channel_list)
mbX1.set_current_ranges(CurrentRange.Range_2mA_SMU, channel_list)
```

```
True
```

We define a simple helper function that can perform a software sweep.

It returns the force values, the measured values and the current range for each step:

```
def sw_sweep(start, stop, step):
    force_values = np.arange(start, stop+step, step)
    ch1_results = np.zeros(len(force_values))
    ch2_results = np.zeros(len(force_values))
    ch1_ranges = []
    ch2_ranges = []
    for i, voltage in enumerate(force_values):
        mbX1.set_voltages(voltage, channel_list)
        measresults = mbX1.measure_channels(wait_for_result=False, sample_count=1, \
↳ repetitions=1, channel_names=channel_list)
        measresult0 = mbX1.get_measurement_results_for_channel(channel1_id)
        ch1_ranges.append(mbX1.idsmu2Modules['M1.S1'].smu.channels[channel1_id].
```

(continues on next page)

(continued from previous page)

```

↪current_range.name)
    ch2_ranges.append(mbX1.idSmu2Modules['M1.S1'].smu.channels[channel2_id].
↪current_range.name)
    ch1_results[i] = measresult0[channel1_id][0]
    ch2_results[i] = measresult0[channel2_id][0]
    return (force_values, ch1_results, ch2_results, ch1_ranges,ch2_ranges)

```

Autoranging in manual IV sweeps

The first step is to perform the sweep with autoranging deactivated. The respective current range is then output for each step.

```

mbX1.set_current_ranges(CurrentRange.Range_200uA_SMU, channel_list)
mbX1.enable_autorange(False, channel_list)
mbX1.set_voltages(1, channel_list)
sweep_without_autorange = sw_sweep(1,4, 0.2)
print (sweep_without_autorange[3])

```

```

['Range_200uA', 'Range_200uA', 'Range_200uA', 'Range_200uA', 'Range_200uA', 'Range_
↪200uA', 'Range_200uA', 'Range_200uA', 'Range_200uA', 'Range_200uA', 'Range_200uA
↪', 'Range_200uA', 'Range_200uA', 'Range_200uA', 'Range_200uA', 'Range_200uA']

```

Autoranging is performed for both channels using the `enable_autorange()` method:

```

mbX1.enable_autorange(True, channel_list)
#Alternatively, autoranging can also be queried or set via the properties/setter of a
↪channel
print (mbX1.idSmu2Modules['M1.S1'].smu.channels[channel1_id].autorange)
sweep_with_autorange = sw_sweep(1,4, 0.2)
print (sweep_with_autorange[4])

```

```

True
['Range_5uA', 'Range_5uA', 'Range_5uA', 'Range_5uA', 'Range_5uA', 'Range_5uA',
↪'Range_5uA', 'Range_20uA', 'Range_200uA', 'Range_200uA', 'Range_200uA', 'Range_
↪200uA', 'Range_2mA', 'Range_2mA', 'Range_2mA', 'Range_2mA']

```

Plotting the results

The following plots show the difference between the sweep with and without autorange. With a fixed value of 200uA, the current curve saturates at approximately this value.

With an autorange set, the better resolution at low currents becomes clear. The curve does not saturate and reaches around 500uA. The plots are annotated with the current ranges:

(Note: the enum prefixes 'Range_' are removed for better readability)

```

from plotly.subplots import make_subplots
def create_fig():
    fig = make_subplots(rows=2, cols=1, subplot_titles=("LED current vs voltage",
↪"LED current vs voltage"))
    fig.add_trace(go.Scatter(x=sweep_without_ourorange[0], y=sweep_without_
↪ourorange[1],
                                mode='lines+markers+text',
                                text = [s.replace("Range_", "") for s in sweep_without_

```


(continues on next page)

(continued from previous page)

```

↪autorange[3]],
        textposition="top center",
        name='LED'), row=1, col=1)
fig.add_trace(go.Scatter(x=sweep_with_outorange[0], y=sweep_with_outorange[1],
        mode='lines+markers+text',
        text = [s.replace("Range_", "") for s in sweep_with_
↪autorange[3]],
        textposition="top center",
        name='LED'), row=2, col=1)
fig.update_yaxes(type="log", range=[np.log(0.0001),np.log(1)], row=1, col=1)
fig.update_yaxes(type="log", range=[np.log(0.0001),np.log(1)], row=2, col=1)
fig.update_layout(height=600)
return fig
fig = create_fig()
fig
#uncomment in pure python script:
#fig.show()

```

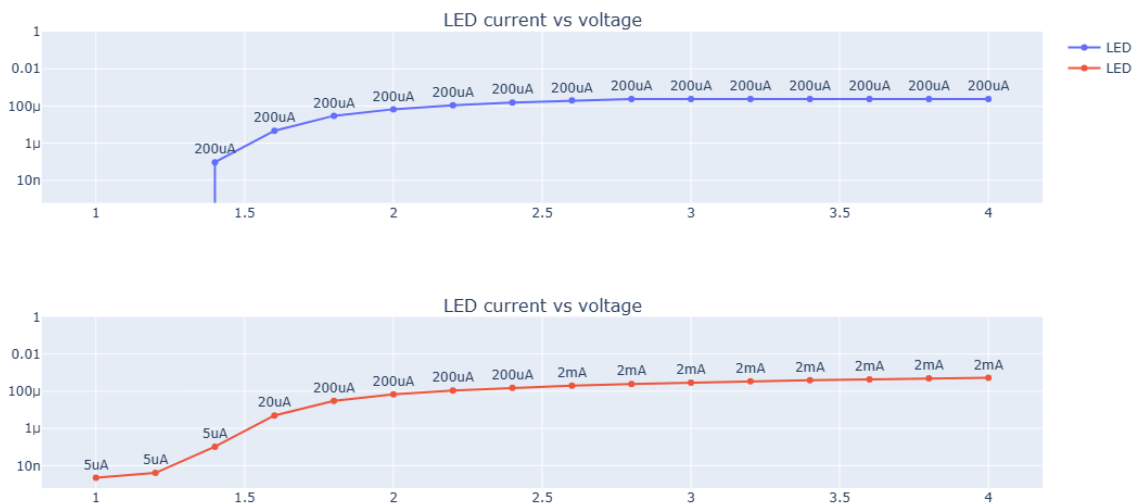

```
srunner.shutdown()
```

0.6 Getting started with the Aspect Device Engine Python API 6

0.6.1 IV Measurements with List Sweeps

This is the sixth introductory overview of programming the Aspect Device Engine Python API.

This document is available as pdf and interactive jupyter notebook. The introduction includes the following objectives:

- Learn how to use the ListSweep class to perform multi-channel sweeps in sub-milliseconds resolution

The *function generator* methods already discussed can be used to generate simple sweeps in real time. With the so-called *list sweeps*, more complex measurements can be carried out simultaneously on several channels of a module.

Please note: In this example, an idSMU board equipped with LEDs with series resistors at the outputs was used. If a differently configured board is used, the measurement curves may look different. In the case of open outputs, the sweep can also be performed via the voltage instead of the current. The voltage is measured back correctly, even with open outputs!

```
from aspectdeviceengine.enginecore import IdSmuService, IdSmuServiceRunner, \
↳ IdSmuBoardModel, MeasurementMode, FunctionGeneratorType, generate_function_
↳ generator_data, CurrentRange
from aspectdeviceengine.enginecore import IdSmuBoardModel, \
↳ ListSweepChannelConfiguration, ListSweep
import plotly.graph_objects as go
from plotly.subplots import make_subplots
import numpy as np
srunner = IdSmuServiceRunner()
mbX1 : IdSmuBoardModel = srunner.get_idsmu_service().get_first_board()
```

```
# It is always a good idea to check if the modules on the board are initialized
print(mbX1.is_board_initialized())
print(mbX1.device_information)
```

```
True
+-----+-----+-----+-----+-----+
| DeviceId | IdSmu-Type | Name   | Firmware | Initialized |
+-----+-----+-----+-----+-----+
| M1.S1    | IdSmu2     | M1.S1  | 0x08191f | true        |
+-----+-----+-----+-----+-----+
```

Preparation

We prepare some variables and enable 2 channels of one idSMU module for the sweep with the multi-side methods. Nothing new here:

```
device_id = "M1.S1"
channel1_id = f'{device_id}.C1'
channel2_id = f'{device_id}.C3'

# assign meaningful alias names for the channels
mbX1.set_channel_name(channel1_id, "LED1")
mbX1.set_channel_name(channel2_id, "LED2")
mbX1.set_device_name(device_id, "My_LED_Module")
channel_list = ['LED1', 'LED2']

# enable the two channels together
mbX1.set_enable_channels(True, channel_list)
idSmu2 = mbX1.idSmu2Modules['M1.S1']

# print some channel information
for i, channel in enumerate(idSmu2.smu.channels.as_list()):
    print(f"Channel {i+1} with name {channel.name}" +
          " is enabled: " + ("YES" if channel.enabled else "NO"))
```

```
Channel 1 with name LED1 is enabled: YES
Channel 2 with name M1.S1.C2 is enabled: NO
Channel 3 with name LED2 is enabled: YES
Channel 4 with name M1.S1.C4 is enabled: NO
```

Preparing the parameters before a sweep

It is **important to understand** that the sweep works with the previously set configurations for the **force type** (voltage for DSP and SMU or current for SMU types) and the **measurement type** (voltage or current). These configurations are *not* part of the parameterisation of the sweep.

With a SMU type there is the option of forcing voltage or current, with a DPS type only voltage. The `set_voltages()` method encapsulates switching to the force type voltage and setting the output value. Although there are also dedicated methods for switching the force type, this is the simplest (for the SMU there is a corresponding `set_currents()` method). To measure a current in the sweep, we set the measurement type to current (*isense*).

We then adjust the current range to the maximum expected current. An autorange does not currently exist for the list sweeps.

```
mbX1.set_voltages(1, channel_list)
mbX1.set_measurement_modes(MeasurementMode.isense, channel_list)
mbX1.set_current_ranges(CurrentRange.Range_2mA_SMU, channel_list) # Note: the ranges_
↪for DPS-Modules differ from SMU modules!
mbX1.get_output_force_value(channel_list[0])
```

```
1.0
```

Configuring the sweep

The `ListSweepChannelConfiguration` and `ListSweep` classes are involved in the parameterisation of a sweep.

`ListSweepChannelConfiguration` is used to configure the form of the sweep for a channel.

`ListSweep` executes the sweep and saves the measurement results for each sweep.

In addition, parameters common to all sweeps are set here, such as the measurement delay.

Single channel sweep

Configuration

We configure a simple linear sweep from 1V to 3V with 20 steps for a single channel. The start and end values are always included in the sweep.

This is why you can expect *number of steps + 1* measurement results.

Important note:

The number of sweep steps is currently limited to 52 steps for a single channel sweep as the memory from which the sweep is executed is limited (see advanced topics)

With 4 lines of code you can configure an executable sweep. In addition, the measurement delay is adjusted in this example:

```
# Instantiation of a channel configuration
config_ch1 : ListSweepChannelConfiguration = ListSweepChannelConfiguration()
```

(continues on next page)

(continued from previous page)

```

# Configuration of a 20-steps linear sweep from 1 to 3 including 1 and 3, in units of
↪ volt
# (because the output force typ was set to voltage before)
config_ch1.set_linear_sweep(1, 3, 20)

# Instantiation of a ListSweep object
sweep_1ch : ListSweep = ListSweep("My_LED_Module", mbX1)

# By adding at least one channel configuration to the list sweep, the sweep is ready
↪ to start.
sweep_1ch.add_channel_configuration("LED1", config_ch1)

# The default measurement delay of 0 (plus some overhead in the lower us-range)
# is often too short. To allow the output voltage to settle before the measurement
# this time is increased to 100 microseconds
sweep_1ch.set_measurement_delay(100)

# let's examine the sweep array before execution:
print(f'Size of sweep array: {len(config_ch1.force_values)}')
print(f'Type of sweep array: {type(config_ch1.force_values)}')
print(f'Values: {config_ch1.force_values}')

```

```

Size of sweep array: 21
Type of sweep array: <class 'numpy.ndarray'>
Values: [1.  1.1 1.2 1.3 1.4 1.5 1.6 1.7 1.8 1.9 2.  2.1 2.2 2.3 2.4 2.5 2.6 2.7
 2.8 2.9 3. ]

```

Execution and results

The sweep is executed by calling the `run()` method.

```

# The run() method starts the sweep and returns after the measurement sweep is
↪ finished
sweep_1ch.run()
# The get_measurement_result() method returns the measured values
currents_LED1 = sweep_1ch.get_measurement_result("LED1")

# The get_force_values() of the configuration object returns the forced values
voltages_LED1 = config_ch1.force_values


# The measurement times (=sample times in case of the sample count = 1 which is
↪ default)
# can be obtained via the sweep object
sample_times = sweep_1ch.timecode

```

Plotting the results

The results are visualised as plots below:

```
fig = make_subplots(rows=1, cols=2, subplot_titles=("LED current vs sample time",
↪ "LED current vs voltage"))
fig.add_trace(go.Scatter(x=sample_times, y=currents_LED1,
                        mode='lines+markers',
                        name='LED1 (Ch1)'), row=1, col=1)
fig.add_trace(go.Scatter(x=voltages_LED1, y=currents_LED1,
                        mode='lines+markers',
                        name='LED1 (Ch1)'), row=1, col=2)
fig.update_xaxes(title_text="Time [us]", row=1, col=1)
fig.update_yaxes(title_text="Current [A]", row=1, col=1)
fig.update_xaxes(title_text="Voltage [V]", row=1, col=2)
fig.update_yaxes(title_text="Current [A]", row=1, col=2)
fig
# uncomment in pure python script:
# fig.show()
```

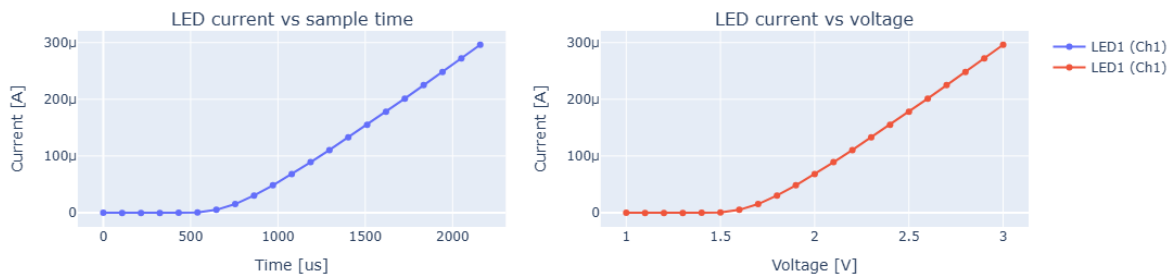

Multi channel sweep

With a few simple modifications for a second channel, a multi-channel sweep is parameterised. The sweeps can be of different lengths. Here we parameterise both channels differently to show the possibility of flexibly setting the sweep parameters for each channel. Instead of the built-in method for a linear sweep, this time we use an array with its own force values on the first channel (even if the numpy arange() function again generates linear values with equidistant values)

```
mbX1.set_current_ranges(CurrentRange.Range_2mA_SMU, channel_list)
# first channel configuration
config_ch1 : ListSweepChannelConfiguration = ListSweepChannelConfiguration()
# custom force values
force_values = np.arange(1, 4, 0.20)
# We modify the numpy array a little to distinguish the sweep curve from the linear_
↪ case
force_values[0:3] = 3
config_ch1.force_values = force_values

# second channel configuration with linear sweep
config_ch2 : ListSweepChannelConfiguration = ListSweepChannelConfiguration()
config_ch2.set_linear_sweep(1, 4, 40)
```

(continues on next page)

(continued from previous page)

```

sweep : ListSweep = ListSweep("My_LED_Module", mbX1)
sweep.add_channel_configuration("LED1", config_ch1)
sweep.add_channel_configuration("LED2", config_ch2)
sweep.set_measurement_delay(100)
sweep.run()

currents_LED1 = sweep.get_measurement_result("LED1")
currents_LED2 = sweep.get_measurement_result("LED2")
sample_times = sweep.timecode

```


Plotting the results

The results are visualised as plots below:

```

fig = make_subplots(rows=1, cols=2, subplot_titles=("LED current vs sample time",
↵ "LED current vs voltage"))
fig.add_trace(go.Scatter(x=sample_times, y=currents_LED1,
                        mode='lines+markers',
                        name='LED1 (Ch1)'), row=1, col=1)
fig.add_trace(go.Scatter(x=sample_times, y=currents_LED2,
                        mode='lines+markers',
                        name='LED2 (Ch3)'), row=1, col=1)
fig.add_trace(go.Scatter(x=config_ch1.force_values, y=currents_LED1,
                        mode='lines+markers',
                        name='LED1 (Ch1)'), row=1, col=2)
fig.add_trace(go.Scatter(x=config_ch2.force_values, y=currents_LED2,
                        mode='lines+markers',
                        name='LED2 (Ch3)'), row=1, col=2)
fig.update_xaxes(title_text="Time [us]", row=1, col=1)
fig.update_yaxes(title_text="Current [A]", row=1, col=1)
fig.update_xaxes(title_text="Voltage [V]", row=1, col=2)
fig.update_yaxes(title_text="Current [A]", row=1, col=2)
fig
#uncomment in pure python script:
#fig.show()

```

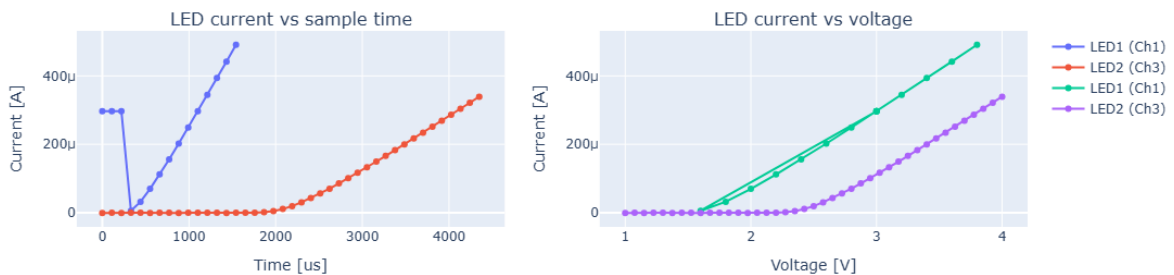

Advanced topics


Changing the current range during a sweep

The list sweep does not support autoranging. However, it is possible to change the range at user-defined points in the sweep.

The last sweep ran with a range of 2mA. We plot the last sweep again and annotate the step indices to the plot.

The switch-on behaviour becomes visible at steps 20-21 or 10-20 microamperes:

```
def create_fig():
    fig = make_subplots(rows=2, cols=1, subplot_titles=("LED current vs sample time",
    ↪ "LED current vs sample time"))
    fig.add_trace(go.Scatter(x=sweep.timecode, y=sweep.get_measurement_result("LED2"),
        mode='lines+markers+text',
        text = [str(i) for i in range(1, len(sample_times) + 1)],
        textposition="top center",
        name='LED1 (Ch1)'), row=1, col=1)
    fig.add_trace(go.Scatter(x=sweep.timecode, y=sweep.get_measurement_result("LED2"),
        mode='lines+markers+text',
        text = [str(i) for i in range(1, len(sample_times) + 1)],
        textposition="top center",
        name='LED1 (Ch1)'), row=2, col=1)
    fig.update_yaxes(type="log", range=[np.log(0.0001), np.log(1)], row=2, col=1)
    fig.update_layout(height=600)
    return fig
fig = create_fig()
fig
#uncomment in pure python script:
#fig.show()
```

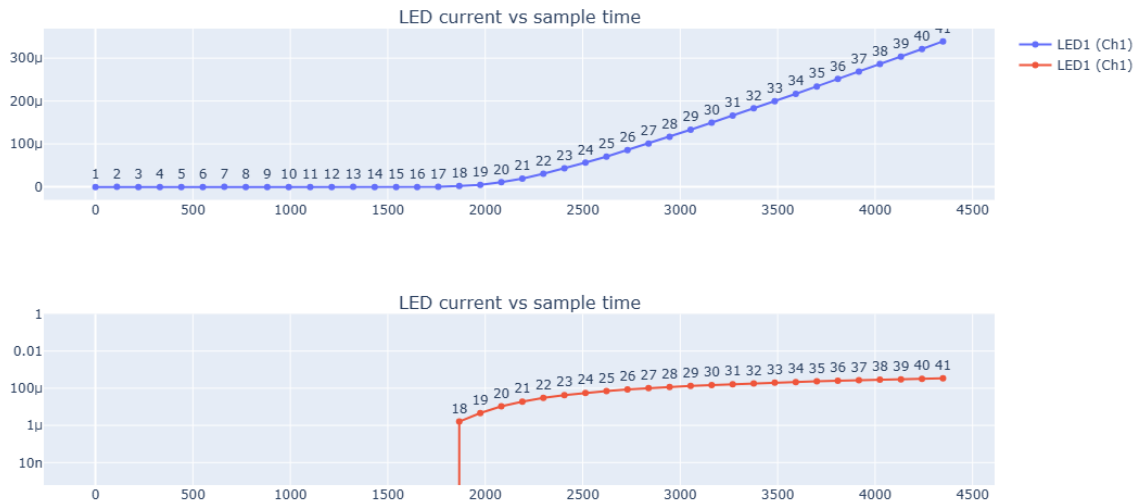

Once you have identified useful points for switching,


you can use the `change_current_range_at()` method to switch the range at these points.

First, we remove any current range switching from the configuration by calling `clear_current_ranges()`.

In this example the sweep starts with a range u 5uA, switches to 20uA at index 18uA and to 200uA at index 20.

As can be seen in the plots below, the range below 20 microamperes is now much better resolved than in the last sweep.

```
mbX1.set_voltages(1, channel_list)
mbX1.set_current_ranges(CurrentRange.Range_5uA, channel_list)
config_ch2.clear_current_ranges()
sweep.set_sample_count(1)
sweep.set_measurement_delay(500)
config_ch2.change_current_range_at(18, CurrentRange.Range_20uA_SMU)
config_ch2.change_current_range_at(20, CurrentRange.Range_200uA_SMU)
sweep.run()
fig = create_fig()
fig
#uncomment in pure python script:
#fig.show()
```

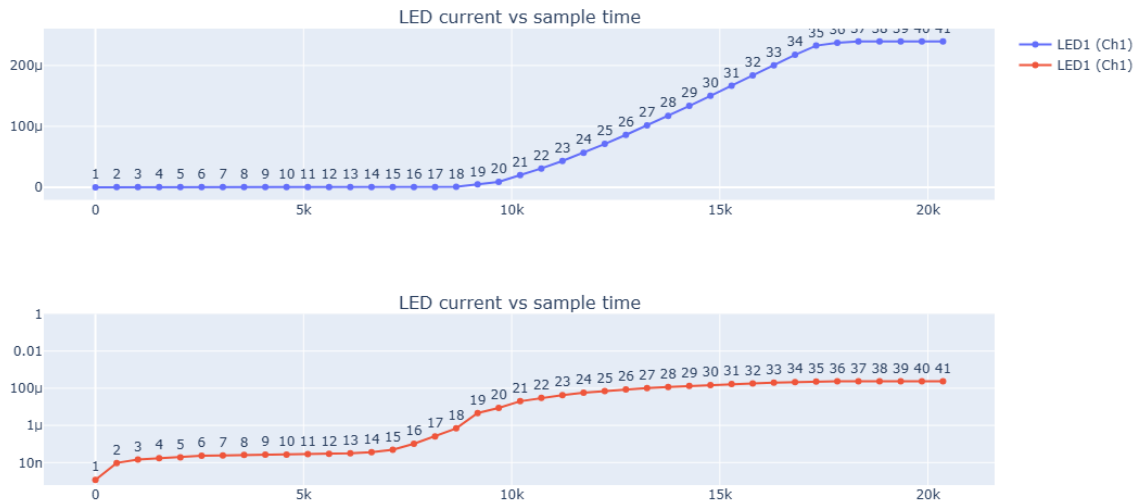

Constant force mode and sweep reusage

Sweeps can be repeated as often as required.

The parameters can remain the same or be changed between two runs.

This is illustrated by the following example where the last sweep is reused.

The first channel is set to 'constant force' mode.

In this mode, only measurements are made and the force value at the output is not varied.

```
config_ch1.set_constant_force_mode(number_of_steps=41)
sweep.run()
currents_LED1 = sweep.get_measurement_result("LED1")
currents_LED2 = sweep.get_measurement_result("LED2")
sample_times = sweep.timecode
```

```
fig = make_subplots(rows=1, cols=1, subplot_titles=("LED current vs sample time",
    "LED current vs voltage"))
fig.add_trace(go.Scatter(x=sample_times, y=currents_LED1,
```

(continues on next page)


(continued from previous page)

```

mode='lines+markers',
name='LED1 (Ch1)'), row=1, col=1)
fig.add_trace(go.Scatter(x=sample_times, y=currents_LED2,
mode='lines+markers',
name='LED2 (Ch3)'), row=1, col=1)

fig.update_xaxes(title_text="Time [us]", row=1, col=1)
fig.update_yaxes(title_text="Current [A]", row=1, col=1)
fig
#uncomment in pure python script:
#fig.show()

```

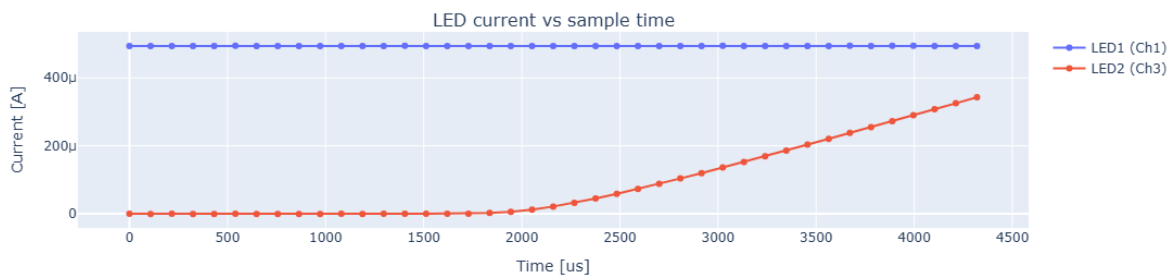

```

runner.shutdown()

```

0.7 Getting started with the Aspect Device Engine Python API 7

0.7.1 Triggering

This is the 7th introductory overview of programming the Aspect Device Engine Python API.

This document is available as pdf and interactive jupyter notebook. The introduction includes the following objectives:

TODO

```

runner.shutdown()

```

0.8 Getting started with the Aspect Device Engine Python API 8a

0.8.1 Parameter tables I

This is the 8th introductory overview of programming the Aspect Device Engine Python API.

This document is available as pdf and interactive jupyter notebook. The introduction includes the following objectives:

- How to load parameter setting tables
- Inspecting the contents of a table
- How to modify the parameters in the table

- Sending the parameters of a table to the hardware

Introduction

There are two basic options for programming the idSMU hardware:

Firstly, by calling object methods or setting properties of the hardware models, as described in the previous chapters.

The second method is a tabular method.

Here, tables are used to hold the values of the hardware parameters. These can be modified and applied.

The tabular method is slower than the direct method calls. However, if performance is not the highest priority, this method can be used to export the state of the hardware and apply it again later.

Instead of dozens of method calls, just one is then sufficient to set all parameters at once.

```
from aspectdeviceengine.enginecore import IdSmuService, IdSmuSettingsService, \
    IdSmuServiceRunner, IdSmuBoardModel, IdSmuSettingsService
from aspectdeviceengine.enginecore import IdSmuBoardModel, IdqTable, IdqTableGroup
import plotly.graph_objects as go
import pathlib, os
import numpy as np
srunner = IdSmuServiceRunner()
mbX1 : IdSmuBoardModel = srunner.get_idsmu_service().get_first_board()
```

The IdSmuSettingsService class is used to load and apply tables:

```
setting_service : IdSmuSettingsService = srunner.get_idsmu_service().get_settings_
    service()
```

Import a table from csv

The parameter tables, also known as setting tables, have a special format.

If you open one of the CSV files with such a table, it looks like this, for example:

[M1_test]												
#SMU-Channel												
SettingId	Hardware	Name	Group	Type	OutputFor	ForceMod	EnableOu	CurrentRa	Measuren	EnableCla	ClampLow	ClampHigh
0	M1.S1.C1	M1.S1.C1	All	SMU-Char	1 FV		1	70mA	HighZ	1	-0,07	0,07
0	M1.S1.C2	M1.S1.C2	All	SMU-Char	1 FV		1	70mA	HighZ	1	-0,07	0,07
0	M1.S1.C3	M1.S1.C3	All	SMU-Char	1 FV		0	70mA	HighZ	1	-0,07	0,07
0	M1.S1.C4	M1.S1.C4	All	SMU-Char	1 FV		0	70mA	HighZ	1	-0,07	0,07
[Types]												
int	string	string	string	string	float	[FV,FI,Hig	bool	[5uA,20uA	[ISense,V'	bool	float	float
#SMU-Device												
SettingId	Hardware	Name	Group	Type	INT10K	ExternalPower						
0	M1.S1	M1.S1	All	SMU-Devi	1	0						
[Types]												
int	string	string	string	string	bool	bool						

Firstly, the name of the table is given in square brackets [].

A table can consist of several groups. As the parameters for the idSMU device are different to the channel parameters, the parameters are each in their own group with their own column headers. The table name is followed by the group name (#).

This is followed by the column header of the group and then the actual data. This is followed by the entries for the data type or, in the case of enumeration types, the valid values for this parameter.

These entries are automatically appended when exporting a table - the user does not have to worry about them.

A table can be loaded with the `import_settings_from_csv()` method of the `SettingService`:

```
setting_service.import_settings_from_csv(os.path.join(os.path.abspath(""), "setting_
←tables_m1.csv" ))
# List all parameter setting names
print(setting_service.get_parameter_settings_names())
# We only print 10 columns so that the output fits on the screen
print(setting_service.print_settings('M1_test', False, False, 10))
```

```

['M1_test']
Group name: SMU-Channel
+-----+-----+-----+-----+-----+-----+-----+
↩+-----+-----+-----+-----+
| SettingId | HardwareId | Name      | Group | Type      | OutputForceValue | FV |
↩+-----+-----+-----+-----+
| 0         | M1.S1.C1   | M1.S1.C1 | All   | SMU-Channel | 1               | FV |
↩      | 1         | 70mA     | HighZ |             |                 |    |
| 0         | M1.S1.C2   | M1.S1.C2 | All   | SMU-Channel | 1               | FV |
↩      | 1         | 70mA     | HighZ |             |                 |    |
| 0         | M1.S1.C3   | M1.S1.C3 | All   | SMU-Channel | 1               | FV |
↩      | 0         | 70mA     | HighZ |             |                 |    |
| 0         | M1.S1.C4   | M1.S1.C4 | All   | SMU-Channel | 1               | FV |
↩      | 0         | 70mA     | HighZ |             |                 |    |
+-----+-----+-----+-----+-----+-----+-----+
↩+-----+-----+-----+-----+
Group name: SMU-Device
+-----+-----+-----+-----+-----+-----+-----+
| SettingId | HardwareId | Name      | Group | Type      | INT10K | ExternalPower |
+-----+-----+-----+-----+-----+-----+-----+
| 0         | M1.S1      | M1.S1     | All   | SMU-Device | 1       | 0              |
+-----+-----+-----+-----+-----+-----+-----+

```

Modification of table entries

The entire table with its groups is managed by an actual table of the type `IdqTable`.

The `dqTableGroup` type represents a group.

```
paratable : IdqTable = setting_service.get_parameter_setting('M1_test')
group : IdqTableGroup = paratable.get_table_group('SMU-Channel')
```

To change a table entry, we need the row index and the column name of the cell.

The helper method `get_row_index()` is used to find the row index for a known entry in the table:

```
row_idx = [group.get_row_index('HardwareId', 'M1.S1.C1'), group.get_row_index(
    ↪ 'HardwareId', 'M1.S1.C2')]
print(row_idx)
```

 $[0, 1]$

The `set_parameter_value()` method is now used to change a value in a table group.

All values in the table are of type string:

```
group.set_parameter_value(row_index=row_idx[0], parameter_name='OutputForceValue',  
↪ parameter_value='5')
```

(continues on next page)

(continued from previous page)

```
group.set_parameter_value(row_index=row_idx[1], parameter_name='OutputForceValue',
↪parameter_value='6')
print(setting_service.print_settings('M1_test', False, False, 10))
```

Group name: SMU-Channel

SettingId	HardwareId	Name	Group	Type	OutputForceValue	
ForceMode	EnableOutput	CurrentRange	MeasurementMode			
0	1	M1.S1.C1	All	SMU-Channel	5	FV
		70mA	HighZ			
0	1	M1.S1.C2	All	SMU-Channel	6	FV
		70mA	HighZ			
0	0	M1.S1.C3	All	SMU-Channel	1	FV
		70mA	HighZ			
0	0	M1.S1.C4	All	SMU-Channel	1	FV
		70mA	HighZ			

Group name: SMU-Device

SettingId	HardwareId	Name	Group	Type	INT10K	ExternalPower
0	M1.S1	M1.S1	All	SMU-Device	1	0

Writing the parameters

The `apply_parameter_setting()` method can now be used to write the parameters of a group or the entire table to the hardware.

We then examine a few channel parameters:

```
setting_service.apply_parameter_setting(setting_name='M1_test', board_address='M1',
↪filtered=False, table_group_name='SMU-Channel')
print(mbX1.idSmu2Modules['M1.S1'].smu.channels['M1.S1.C1'].enabled)
print(mbX1.idSmu2Modules['M1.S1'].smu.channels['M1.S1.C1'].voltage)
print(mbX1.idSmu2Modules['M1.S1'].smu.channels['M1.S1.C2'].enabled)
print(mbX1.idSmu2Modules['M1.S1'].smu.channels['M1.S1.C2'].voltage)
```

```
True
4.998809814453125
True
5.999847412109375
```

```
group.at[row_idx[0], 'OutputForceValue'] = "0"
group.at[row_idx[1], 'EnableOutput'] = "0"
print(setting_service.print_settings('M1_test', False, False, 10))
setting_service.apply_parameter_setting('M1_test', 'M1', False, 'SMU-Channel')
```

```

Group name: SMU-Channel
+-----+-----+-----+-----+-----+-----+-----+
↪ | SettingId | HardwareId | Name      | Group | Type      | OutputForceValue |
↪ | ForceMode | EnableOutput | CurrentRange | MeasurementMode |
+-----+-----+-----+-----+-----+-----+-----+
↪ | 0         | M1.S1.C1   | M1.S1.C1 | All   | SMU-Channel | 0                | FV ↪
↪ | 1         |             | 70mA     | HighZ |             |                  |
↪ | 0         | M1.S1.C2   | M1.S1.C2 | All   | SMU-Channel | 6                | FV ↪
↪ | 0         |             | 70mA     | HighZ |             |                  |
↪ | 0         | M1.S1.C3   | M1.S1.C3 | All   | SMU-Channel | 1                | FV ↪
↪ | 0         |             | 70mA     | HighZ |             |                  |
↪ | 0         | M1.S1.C4   | M1.S1.C4 | All   | SMU-Channel | 1                | FV ↪
↪ | 0         |             | 70mA     | HighZ |             |                  |
+-----+-----+-----+-----+-----+-----+-----+
↪ +-----+-----+-----+-----+-----+-----+-----+

Group name: SMU-Device
+-----+-----+-----+-----+-----+-----+-----+
| SettingId | HardwareId | Name  | Group | Type      | INT10K | ExternalPower |
+-----+-----+-----+-----+-----+-----+-----+
| 0         | M1.S1      | M1.S1 | All   | SMU-Device | 1      | 0              |
+-----+-----+-----+-----+-----+-----+-----+

```

4

```
srunner.shutdown()
```

0.9 Getting startet with the Aspect Device Engine Python API 8b

0.9.1 Parameter tables II

This is the 8th introductory overview of programming the Aspect Device Engine Python API.

This document is available as pdf and interactive jupyter notebook. The introduction includes the following objectives:

- How to export the state of the hardware into a parameter table file
- Filtering

Introduction

The basics of the parameter tables were presented in the last chapter.

Further options will now be discussed.

```

from aspectdeviceengine.enginecore import IdSmuService, IdSmuSettingsService, ↪
↪ IdSmuServiceRunner, IdSmuBoardModel, IdSmuSettingsService
from aspectdeviceengine.enginecore import IdSmuBoardModel, IdqTable, IdqTableGroup
import plotly.graph_objects as go
import pathlib, os
import numpy as np
srunner = IdSmuServiceRunner()
mbX1 : IdSmuBoardModel = srunner.get_idsmu_service().get_first_board()

```

The `IdSmuSettingsService` class is used to load and apply tables:

```
setting_service : IdSmuSettingsService = srunner.get_idsmu_service().get_settings_
    ↵service()
```

Exporting a table

The Methode `get_parameter_settings_for_board()` is used to read out the current status of the hardware and generate a parameter table.

We first change a few parameters and then analyse the table:

```
mbX1.idSmu2Modules['M1.S1'].smu.channels['M1.S1.C1'].voltage = 1
mbX1.idSmu2Modules['M1.S1'].smu.channels['M1.S1.C2'].voltage = 2
mbX1.idSmu2Modules['M1.S1'].smu.channels['M1.S1.C3'].voltage = 3
mbX1.idSmu2Modules['M1.S1'].smu.channels['M1.S1.C4'].voltage = 4
mbX1_settings : IdqTable = setting_service.get_parameter_settings_for_board('M1')
print(mbX1_settings.name)
print(setting_service.print_settings(mbX1_settings.name, False, False, 10))
```

```

M1
Group name: SMU-Channel
+-----+-----+-----+-----+-----+-----+-----+
↩+-----+-----+-----+-----+
| SettingId | HardwareId | Name      | Group | Type      | OutputForceValue | FV |
↩+-----+-----+-----+-----+-----+
| ForceMode | EnableOutput | CurrentRange | MeasurementMode |
+-----+-----+-----+-----+-----+
↩+-----+-----+-----+-----+
| 0          | M1.S1.C1    | M1.S1.C1 | All   | SMU-Channel | 1,000000         | FV |
↩          | 0           | 70mA     | HighZ |              |                   |    |
| 0          | M1.S1.C2    | M1.S1.C2 | All   | SMU-Channel | 2,000000         | FV |
↩          | 0           | 70mA     | HighZ |              |                   |    |
| 0          | M1.S1.C3    | M1.S1.C3 | All   | SMU-Channel | 3,000000         | FV |
↩          | 0           | 70mA     | HighZ |              |                   |    |
| 0          | M1.S1.C4    | M1.S1.C4 | All   | SMU-Channel | 4,000000         | FV |
↩          | 0           | 70mA     | HighZ |              |                   |    |
+-----+-----+-----+-----+-----+
↩+-----+-----+-----+-----+

Group name: SMU-Device
+-----+-----+-----+-----+-----+-----+-----+
| SettingId | HardwareId | Name      | Group | Type      | INT10K | ExternalPower |
+-----+-----+-----+-----+-----+-----+-----+
| 0          | M1.S1      | M1.S1     | All   | SMU-Device | 1       | 0              |
+-----+-----+-----+-----+-----+-----+

```

Before saving the table, it is advisable to change the name of the table.

The default name is the board address.

The name of a table to be loaded later should therefore be different from this default value so that this table is not overwritten by the method executed above.

```
mbX1_settings.name = "M1_voltages_set"
print(mbX1_settings.name)
print(setting_service.get_parameter_settings_names())
```

```
M1_voltages_set
['M1_voltages_set']
```

The table is now exported using the `export_settings_to_csv()` method. The three parameters of the method are self-explanatory.

The reason why `setting_names` is a list is that it is also possible to write several tables to one file.

In the example here, we only write the table just created to the file:

```
setting_service.export_settings_to_csv(file_path=os.path.join(os.path.abspath(""),
                                                             mbX1_settings.name + ".csv"), setting_names=[
    ↪ 'M1_voltages_set'],
                                     append_to_file=False)
```

The result looks like this:

[M1_voltages_set]												
#SMU-Channel												
SettingId	Hardware	Name	Group	Type	OutputForceValue	ForceMod	EnableOu	CurrentRa	Measuren	EnableCla	ClampLow	ClampHighValue
0	M1.S1.C1	M1.S1.C1	All	SMU-Char	1	FV		0 70mA	HighZ	1	-0,07	0,07
0	M1.S1.C2	M1.S1.C2	All	SMU-Char	2	FV		0 70mA	HighZ	1	-0,07	0,07
0	M1.S1.C3	M1.S1.C3	All	SMU-Char	3	FV		0 70mA	HighZ	1	-0,07	0,07
0	M1.S1.C4	M1.S1.C4	All	SMU-Char	4	FV		0 70mA	HighZ	1	-0,07	0,07
[Types]												
int	string	string	string	string	float	[FV,FI,Hig bool		[5uA,20uA [ISense,V' bool			float	float
#SMU-Device												
SettingId	Hardware	Name	Group	Type	INT10K	ExternalPower						
0	M1.S1	M1.S1	All	SMU-Devi	1	0						
[Types]												
int	string	string	string	string	bool	bool						

Filtering and applying sub-tables

Previously, all entries in a table or group were always sent to the hardware.

But what if you only want to use certain entries?

There are various filter methods. You can then send only the filtered table to the hardware instead of the entire table.

Row filter

The parameters of the `filter_rows()` method are the column name, the value in this column and whether to search for whole words and not just substrings (`exact_match`).

To print the filtered version of a table, the second parameter in the `print_settings()` method is set to `True`.

```
filtered_table = mbX1_settings.filter_rows(column_name='HardwareId', filter_value='M1.
    ↪ S1.C2', exact_match=False)
print(setting_service.print_settings('M1_voltages_set', True, False, 10))
```

```
Group name: SMU-Channel
+-----+-----+-----+-----+-----+-----+-----+-----+-----+-----+
    ↪ +-----+-----+-----+-----+-----+-----+-----+-----+-----+-----+
    | SettingId | HardwareId | Name      | Group | Type      | OutputForceValue | ForceMod | EnableOutput | CurrentRange | MeasurementMode |
    ↪ +-----+-----+-----+-----+-----+-----+-----+-----+-----+-----+
    | 0         | M1.S1.C2   | M1.S1.C2 | All   | SMU-Channel | 2,000000         | FV       |             | 70mA         | HighZ           |
    ↪ | 0         |             |           |       |             |                   |          |             |              |                 |
    +-----+-----+-----+-----+-----+-----+-----+-----+-----+-----+
    ↪ +-----+-----+-----+-----+-----+-----+-----+-----+-----+-----+
```

The pipe character (|) can be used to realise a logical 'or':


```

filtered_table = mbX1_settings.filter_rows(column_name='HardwareId', filter_value='M1.
↳S1.C2|M1.S1.C4', exact_match=False)
print(setting_service.print_settings('M1_voltages_set', True, False, 10))

```

```

Group name: SMU-Channel
+-----+-----+-----+-----+-----+-----+-----+
↳+-----+-----+-----+-----+-----+
| SettingId | HardwareId | Name      | Group | Type      | OutputForceValue | FV |
↳+-----+-----+-----+-----+-----+
| 0         | M1.S1.C2   | M1.S1.C2 | All   | SMU-Channel | 2,000000         | FV |
↳| 0         |             | 70mA     | HighZ |             |                   |    |
| 0         | M1.S1.C4   | M1.S1.C4 | All   | SMU-Channel | 4,000000         | FV |
↳| 0         |             | 70mA     | HighZ |             |                   |    |
+-----+-----+-----+-----+-----+-----+
↳+-----+-----+-----+-----+-----+

```

If the 'filtered' flag of the 'apply_parameter_setting()' method is set, only the filtered version is written to the hardware instead of the entire table:

```

mbX1.idSmu2Modules['M1.S1'].smu.channels['M1.S1.C1'].voltage = 5
mbX1.idSmu2Modules['M1.S1'].smu.channels['M1.S1.C2'].voltage = 5
mbX1.idSmu2Modules['M1.S1'].smu.channels['M1.S1.C3'].voltage = 5
mbX1.idSmu2Modules['M1.S1'].smu.channels['M1.S1.C4'].voltage = 5
setting_service.apply_parameter_setting(setting_name='M1_voltages_set', board_address=
↳'M1', filtered=True, table_group_name='SMU-Channel')

```

2

We export the status of the hardware to a new table. Now only the values of the filtered table should have been written. The voltages of channels 1 and 3 were set to 5 volts. The remaining channels retain their values.

```

mbX1_settings : IdqTable = setting_service.get_parameter_settings_for_board('M1')
print(setting_service.print_settings(mbX1_settings.name, False, False, 10))

```

```

Group name: SMU-Channel
+-----+-----+-----+-----+-----+-----+-----+
↳+-----+-----+-----+-----+-----+
| SettingId | HardwareId | Name      | Group | Type      | OutputForceValue | FV |
↳+-----+-----+-----+-----+-----+
| 0         | M1.S1.C1   | M1.S1.C1 | All   | SMU-Channel | 5,000000         | FV |
↳| 0         |             | 70mA     | HighZ |             |                   |    |
| 0         | M1.S1.C2   | M1.S1.C2 | All   | SMU-Channel | 2,000000         | FV |
↳| 0         |             | 70mA     | HighZ |             |                   |    |
| 0         | M1.S1.C3   | M1.S1.C3 | All   | SMU-Channel | 5,000000         | FV |
↳| 0         |             | 70mA     | HighZ |             |                   |    |
| 0         | M1.S1.C4   | M1.S1.C4 | All   | SMU-Channel | 4,000000         | FV |
↳| 0         |             | 70mA     | HighZ |             |                   |    |
+-----+-----+-----+-----+-----+-----+
↳+-----+-----+-----+-----+-----+

```

Group name: SMU-Device

```

+-----+-----+-----+-----+-----+-----+-----+

```

(continues on next page)

(continued from previous page)

SettingId	HardwareId	Name	Group	Type	INT10K	ExternalPower
0	M1.S1	M1.S1	All	SMU-Device	1	0

```
srunner.shutdown()
```